


## Leak Test Procedure

### To perform leak test...

- Close APL (Scavenging /Adjustable Pressure Limiting) valve by turning knob clockwise.
- Place thumb over patient connection of breathing circuit Y.
- Remove breathing bag and cover bag port opening. (Use palm of hand that is covering Y.)
- With oxygen (50-55 PSI [3.4-3.8 Bar]) supplied to anesthesia machine, slowly open flowmeter to register 30cm H<sub>2</sub>O on anesthesia machine pressure gauge.
- Turn off flowmeter when pressure reaches 30cm H<sub>2</sub>O. (If pressure holds steady the system is leak free but if pressure drops, proceed to step (F).)
- Slowly open flowmeter until pressure stabilizes at 30cm. H<sub>2</sub>O setting. (This determines the magnitude of the leak. If leak rate is greater than 300ml/min; proceed to step (G).)
- Refer to page 2, "What if machine leaks?"
- Replace reservoir bag. Repeat step (B) and steps (D) through (F). This will determine the integrity of breathing bag.

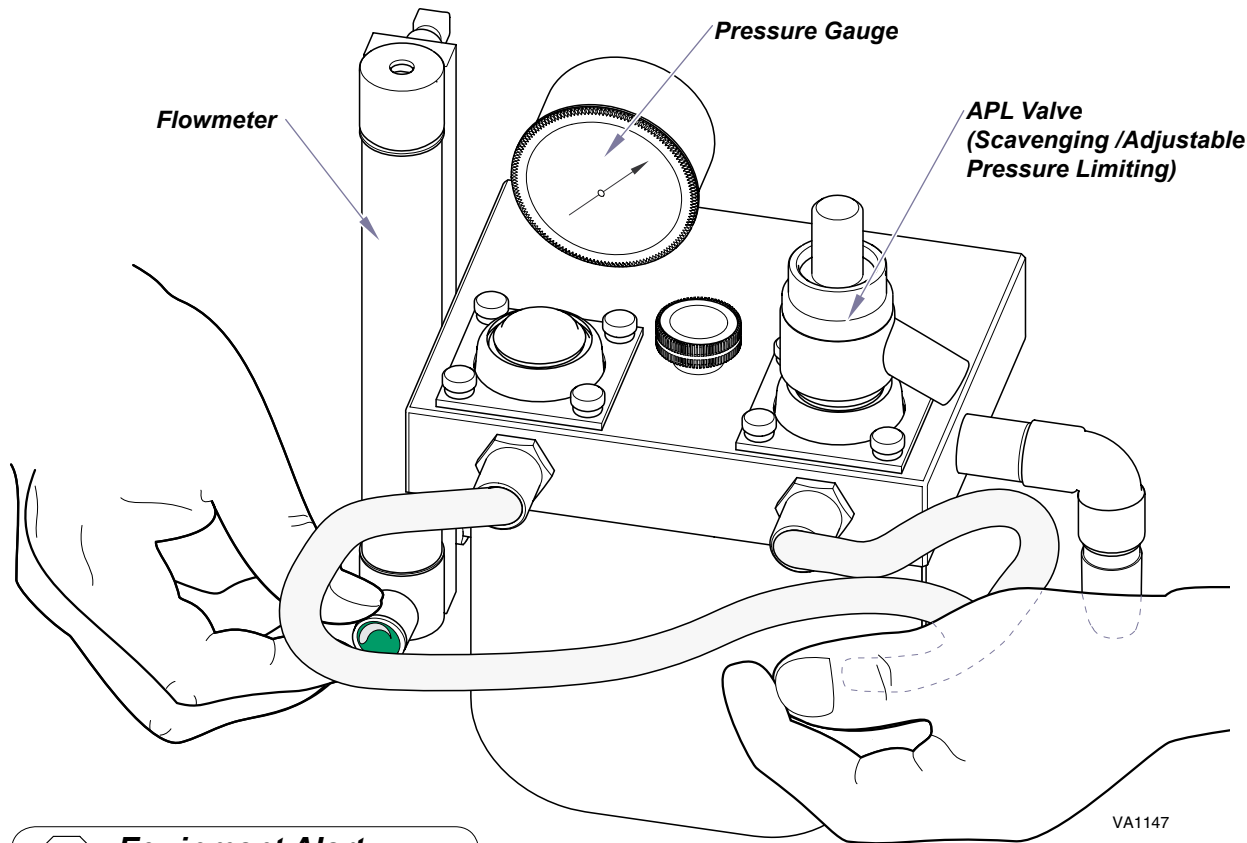


VA1147



### Equipment Alert

Do not activate the oxygen flush during any part of this leak test.

## Maintenance Kit

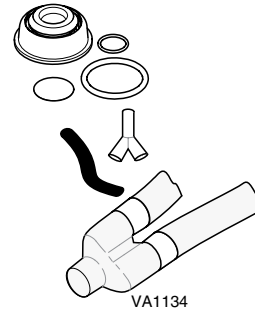

### VME™ Tabletop & VME2®

Order #: 91303017

Includes:

- Disc (2)
- Large and Small O-rings (2 of each size)
- Domes
- 2 Liter Breathing Bag
- Tubing & Y Connector
- Breathing Circuit
- Reservoir Bag
- Gasket

### Sodasorb®

Order #: 002-1761-00

Sodasorb® LF, Canister Pak, Case

Order #: 002-1762-00

Sodasorb® Standard, Canister Pak, Case

Order #: 002-1763-00

Sodasorb® LF, Pre-pak, Case

Order #: 002-1764-00

Sodasorb® Standard, Pre-pak, Case

Order #: 002-1765-00

Sodasorb® Standard, Bucket

### Calling for Service

If service is required, contact your Midmark dealer.

To contact Midmark directly:

8:00am until 5:00 PM Monday through Friday (EST)

**Customer Service**

1-800-Midmark (1-800-643-6275), Fax 1-877-725-6495

**Technical Service**

1-888-279-1260, Fax 1-716-662-8440

### Registration

To register your product warranty, go to [www.midmark.com](http://www.midmark.com)

MIDMARK®

## Quick Reference Guide Matrix VME™ Tabletop & VME2® Anesthesia Machines

Model Number:

Serial Number:

Dealer Information:

Purchase Date:

10578500 Rev C

The Operation of Matrix® anesthesia machines is restricted to use by veterinary professionals trained in anesthesia practice.

This guide outlines basic operation and maintenance of the Anesthesia Machine. Additional information is available on [www.documark.com](http://www.documark.com). If service is required, contact your authorized Midmark dealer.

MIDMARK®




### CAUTION

Insure gas supplies are adequate and turned on. Verify vaporizer is filled to eliminate the possibility of anesthetic depletion during surgery.

#### To Operate Anesthesia Machine...

- A) Loosen thumb screw to remove CO2 absorber canister. Fill canister with absorbent material, following directions on canister label. Be sure canister and gasket mating surface are completely free of absorbent. Replace canister in absorber assembly. Tighten thumb screw.
- B) Connect breathing bag and breathing circuit.
- C) Connect a 50-55 PSI (3.4-3.8 Bar) oxygen supply line to the oxygen supply fitting on the back of the flowmeter.
- D) Adjust flowmeter and vaporizer settings to meet physiologic needs of the patient.
- E) During anesthesia, monitor the pressure gauge, inhalation and exhalation valves, and the breathing bag. Make necessary corrections in flow rate, vaporizer setting, and Scavenging/Adjustable Pressure Limiting (APL) valve, to insure adequate depth of anesthesia and adequate ventilation of the patient.

Note: For guidelines on induction and anesthesia techniques, or particular species physiologic requirements, etc., please refer to a veterinary anesthesia textbook.

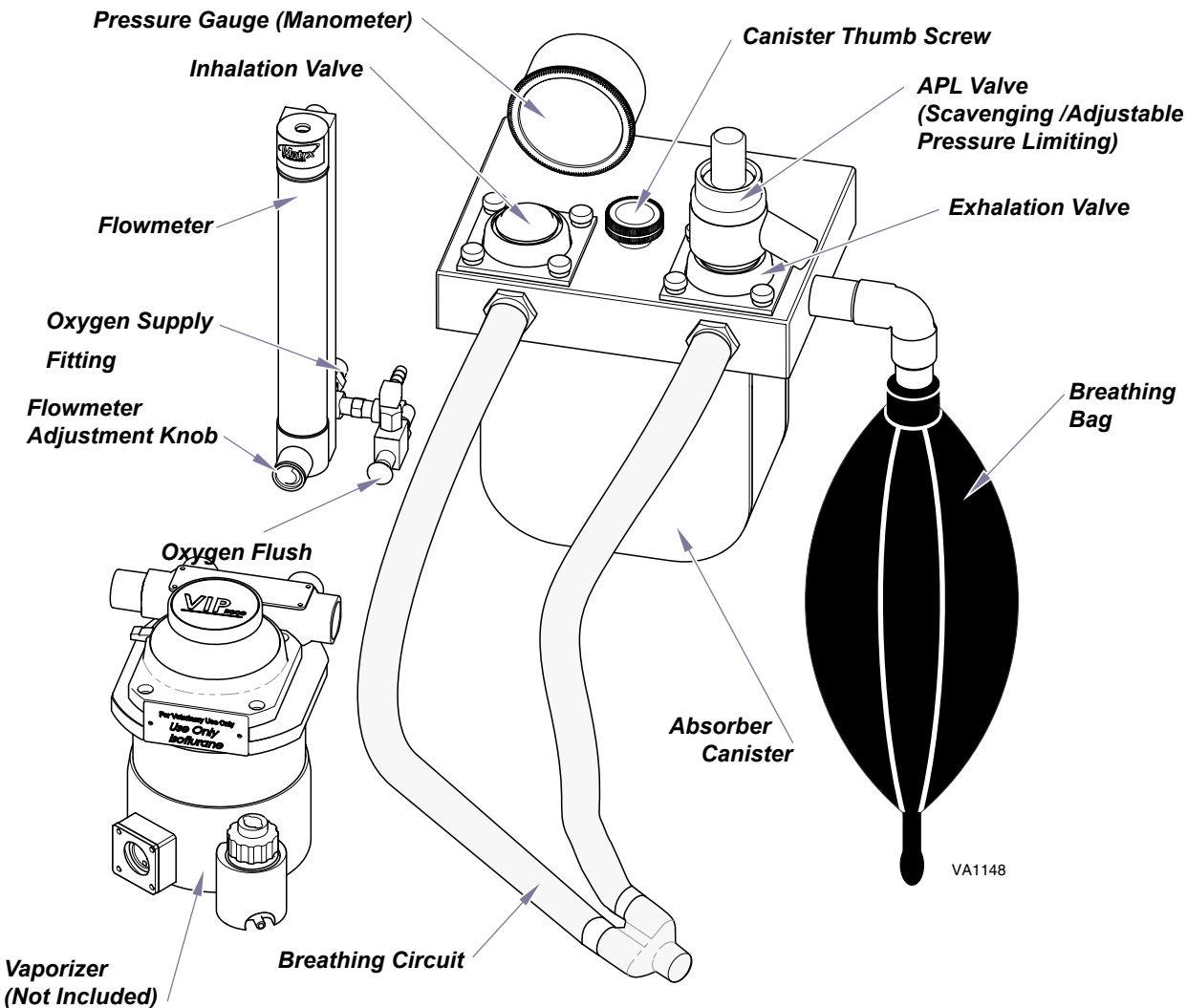

| BREATHING BAG         |          |
|-----------------------|----------|
| Patient Weight        | Bag Size |
| 15 lbs or less (7kg)  | 500 ml   |
| 15-30 lbs (7-14 kg)   | 1 liter  |
| 30-60 lbs (14-27 kg)  | 2 liter  |
| 60-100 lbs (27-45 kg) | 3 liter  |
| 100 + lbs (45 + kg)   | 5 liter  |

#### What if Machine Leaks?

- 1) **Breathing Bag** - If leak occurs, replace.
- 2) **Breathing Circuit** - Install new breathing circuit or obstruct inhalation / exhalation openings to determine if leak originates from breathing circuit.
- 3) **Vaporizer Fittings** - Verify fittings and tubing are securely attached.
- 4) **Canister Gasket** - Check for loose absorbent grains between canister housing gasket.
- 5) **Verify Canister** is seated properly.
- 6) **APL Valve** - Remove valve and obstruct opening to determine if leak originates from APL valve. Check the o-ring under valve for damage.
- 7) **O-Rings** under dome and clamp. Check the o-rings for damage.

## Maintenance Schedule

### Before Every Surgery

#### A. Inspect.

Check machine connections and rubber parts for looseness, damage or wear. Replace as necessary.

#### B. Perform leak test.

### Daily

#### A. Clean.

Remove breathing bag and breathing circuit. Wash with warm water and mild soap, rinse well. Hang to dry. Remove white disks from inhalation and exhalation valves. Wipe with soft cloth, wipe out valves and reassemble. Remove absorber canister, wipe off gaskets, canister and absorber housings.

### Weekly

#### A. Change Absorbent Material.

Follow manufacturer's recommendations for proper use of CO2 absorbents.

### Annually

#### . Service & Calibrate Vaporizer

Contact Dealer for details.

### Every Two (2) Years

#### A. Replace Maintenance Components.

Replace components found in maintenance kit 91303017.



#### Equipment Alert

Never invert or tilt the vaporizer when it contains anesthetic liquid. Never overtighten the flowmeter needle valves, damage to the needle and seat will eventually result.